

NG3000A

générateur d'azote

Numéro de pièce: 10-1210 (230v) - UK |
10-1220 (230v) - EU | 10-1230 (110v) -
US



La description

Le NG3000A fournit un azote d'un ultra haut degré de pureté à des limites de détection de traces grâce à la technologie d'absorption à balancement de pression à des débits de 3 L/min pour les applications GC.



Types d'application


GC, ICP, DSC, DMA

Generator Details

- **Types de gaz:** Azote
- **Compresseur intégré:** Oui

Principales caractéristiques

- Convient aux applications telles que GC, DSC, DMA
- Convient aux applications nécessitant une très haute pureté
- Disponible avec ou sans compresseur externe
- Contient une colonne de tamis moléculaire de carbone auto-régénérant
- Le gaz est fourni sur demande afin que le générateur fonctionne selon votre horaire
- Garantie complète sur site de 12 mois

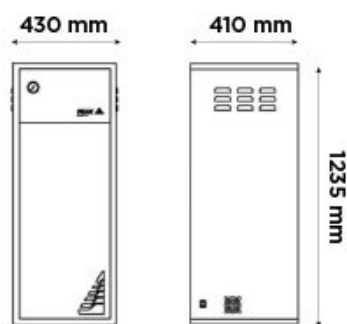

Spécifications techniques

	Générateur d'azote NG3000A
Type de gaz	Azote
Débit de gaz max:	3
Pression de sortie max:	80psi/5,5bar
Pression de sortie max:	80psi/5,5bar
Pureté max:	99.9995%
Voltage:	115 / 230 ± 10% VAC
Fréquence:	50 / 60 Hz
Courant:	3.6 Amps / 8 Amps
Accréditations:	CE
Taille (HxLxP) mm	1235 x 449 x 406,5 mm
Taille (HxLxP) pouces	48,6 x 17,7 x 16 inches
Poids du générateur	88kg / 193,6lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

