

NG3000

générateur d'azote

Numéro de pièce: 10-2210 (230v) - UK |
10-2211 (230v) - EU | 10-2230 (110v) -
US




La description

Le NG3000 fournit un azote d'un ultra haut degré de pureté à des limites de détection de traces grâce à la technologie d'absorption à balancement de pression à des débits de 3 L/min pour les applications GC.



Generator Details

- **Types de gaz:** Azote
- **Compresseur intégré:** Non

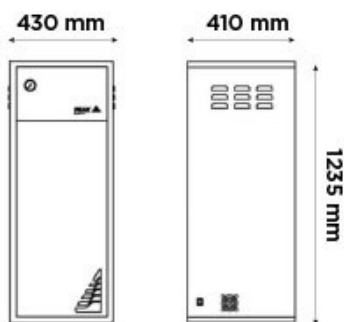

Principales caractéristiques

- Convient aux applications telles que GC, DSC, DMA
- Convient aux applications nécessitant une très haute pureté
- Disponible avec ou sans compresseur externe
- Contient une colonne de tamis moléculaire de carbone auto-régénérant
- Le gaz est fourni sur demande afin que le générateur fonctionne selon votre horaire
- Garantie complète sur site de 12 mois

Spécifications techniques

Générateur d'azote NG3000

Type de gaz	Azote
Débit de gaz max:	3
Pression de sortie max:	80psi/5,5bar
Pression de sortie max:	80psi/5,5bar
Pression d'entrée min:	120psi/8,3bar
Pression d'entrée max:	130psi/9bar
Voltage:	110 / 230 ± 10% VAC
Fréquence:	50 / 60 Hz
Courant:	0.5 Amps / 1.2 Amps
Température de fonctionnement max:	25°C / 77°F
Accréditations:	CE
Taille (HxLxP) mm	1260 x 430 x 410 mm
Taille (HxLxP) pouces	49,6 x 16,9 x 16,1 inches
Poids du générateur	77kg / 169,4lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

